


# get to know FreedomEdge™

## Syringe Infusion System

### Easy to use & learn

No controls - just load, close and go!

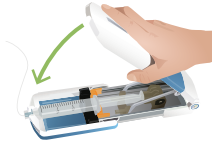

### Silent operation

Convenient for patients - no unnecessary alarms & noise.



### Safe constant pressure

Highly accurate constant pressure pump for safe delivery. Maintains pressure of 13.5 psi.



### Ideal for all users

Maximum simplicity for all users - small hands, limited dexterity and limited strength.



### DynEQ™ (Dynamic Equilibrium)

If the pressure in your sites increases, the system will immediately respond by automatically slowing down.



### precision FLOW RATE TUBING™

RMS Precision Flow Rate Tubing™ controls the infusion rate in the FreedomEdge™ system.




### Portable

No batteries or electricity needed. Comes with a travel pouch to promote patient independence.



### Luer disc

Helps prevent contamination. For patient safety, the pump is designed not to work unless a set with the disc is present.

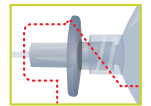


### Quality assured

Made in the USA with the highest level of quality control for defect-free pumps. No product maintenance needed.



### Deliver in minutes to days

Sets from F0.5 to F2400 offer precise infusion rates for a variety of administration times.

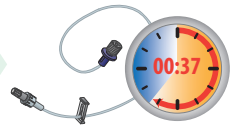


### Fits 20 & 30 ml syringes

Easy to handle and purge of air. Can be filled from 1ml to 30ml.



### Emergency slide clamp

Stop flow immediately, if needed. F-number printed on clamp to assure correct tubing is being used.

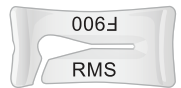

### Progress window

Easy to check infusion progress.



### Low residual volume

Consistent delivery of medication for a full therapeutic dose without expensive drug waste, ≤ 0.16ml.



rmsmedicalproducts.com  
845-469-2042 800-624-9600

The FreedomEdge™ Syringe Infusion System, RMS Precision Flow Rate Tubing™ and DynEQ™ are registered trademarks of RMS Medical Products and are compliant with Medical Device Directive 93/42/EEC. RMS Medical Products is ISO 13485 certified. © 2016 RMS Medical Products; All Rights Reserved. RHI\_INFG\_FEbenefits\_03-10007vA