

get to know HlgH-Flo™

Subcutaneous Safety Needle Sets

Backcut needle tip design

Minimizes damage, scarring & local site reactions. Promotes faster healing.




Extension set

Available for those who prefer longer length.



Flexible flying hinge

Wings have great flexibility & allow each needle to stay in place without any tenting.

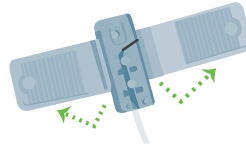


Y-connector

Combines sets for more than 6 sites (up to 8), or to enable infusing with different length needles.



Safety closure

Protects against accidental needle sticks.

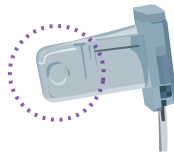


3M Tegaderm™

Transparent dressing that is preferred by most patients is included with all sets.



Custom luer



Designed to provide even flow to every site.

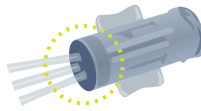

26g HlgH-Flo Needle Sets Options (Residual Volume)


Gauge options

26g needles designed for sub-q flow rates up to 40ml/hr/site & 24g for up to 300ml/hr/site.

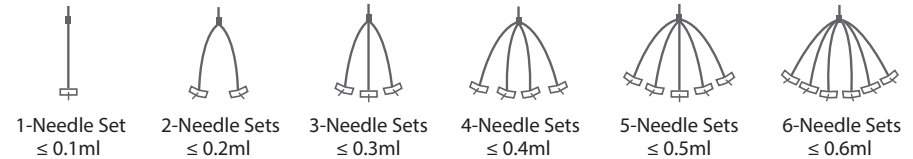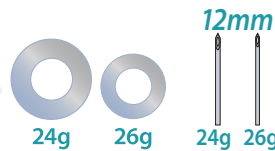

Length options


Needles are available in 4mm, 6mm, 9mm, 12mm, & 14mm lengths.

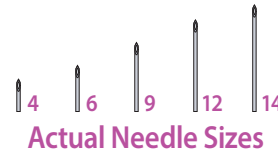

24g HlgH-Flo Needle Sets Options (Residual Volume)

Tubing length

Efficient & practical 20" length (in 26g sets) for less tangles & faster flow rates. 24" in 24g sets.

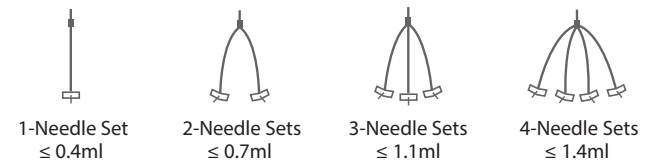

rmsmedicalproducts.com
845-469-2042 800-624-9600