


get to know precision™ Flow Rate Tubing

Part of the Freedom Integrated Syringe Infusion System

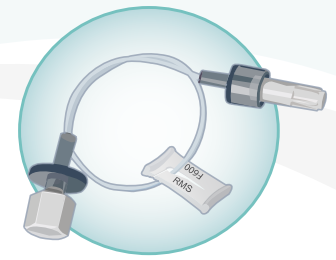


Safe and reliable

Precision is designed for uniform flow profiles, with no free-flow, bolus or overdose.



Accurate infusion rates

Low residual volumes & precise dimensions help to achieve the desired flow rate.

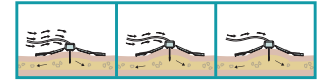

DynEQ™ (Dynamic Equilibrium)

If the pressure in your sites increases, the system will immediately respond by automatically slowing down.



Quality assured

Made in the USA with the highest level of quality control for defect-free products. No maintenance needed.




Low residual volumes

Consistent delivery of medication for a full therapeutic dose without expensive drug waste, $\leq 0.16\text{ml}$.




Emergency slide clamp

Stop flow immediately, if needed. F-number printed on clamp to assure correct tubing is being used.

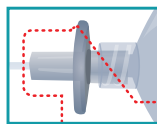

Luer disc

Helps prevent contamination. For patient safety, the syringe driver is designed not to work unless a set with the disc is present.

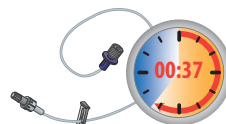

Deliver in minutes to days

Sets from F0.5 to F2400 offer precise infusion rates for a variety of administration times.



RMS flow rate calculator

Convenient reference tool to help determine the appropriate tubing set for desired sub-q infusion duration (rmscalculator.com).

Residual volumes: very low, low and high flow

- **Very low flow** $\leq 0.14 \text{ ml}$: F0.5, F1, F2, F3, F3.8, F5, F8, F10, F15
- **Low flow** $\leq 0.16 \text{ ml}$: F30, F45, F60, F120, F180
- **High flow** $\leq 0.15 \text{ ml}$: F275, F420, F500, F600, F900, F1200, F2400

NOTE: For use with a RMS Freedom Syringe Driver and RMS Freedom System accessories. Use with any other brand may cause over delivery of fluids, unknown flow rates or additional site complications such as pain, swelling and redness.

Related Freedom Integrated System Products:



FREEDOM60® Syringe Infusion System

For larger infusion doses, 60ml syringe



FreedomEdge® Syringe Infusion System

For smaller infusion doses, 20ml & 30ml syringes



HiGH-Flo Subcutaneous Safety Needle Sets™

24G & 26G, 4 -14mm lengths, 1-6 leg sets



24" Extension Tubing Set (FEXT)

For those who need longer tubing length



rmsmedicalproducts.com 845-469-2042 800-624-9600

The FREEDOM60® and FreedomEdge® Syringe Infusion Systems, RMS Precision Flow Rate Tubing™, HiGH-Flo Subcutaneous Safety Needle Sets™, and DynEQ™ are registered trademarks of RMS Medical Products and are compliant with Medical Device Directive 93/42/EEC. RMS Medical Products is ISO 13485 certified. © 2017 RMS Medical Products; All Rights Reserved. RHL_INFG_PRCbenefits_03-10006VC