
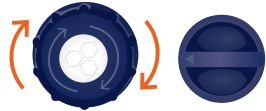



# get to know FREEDOM60<sup>®</sup> Syringe Driver

Part of the Freedom Integrated Syringe Infusion System

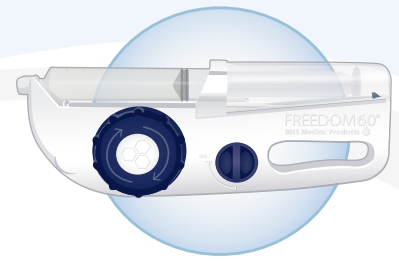


## Easy to use & learn

Only 2 controls.  
Just wind, load & go!

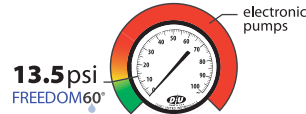

## Safe constant pressure

Highly accurate constant pressure system for safe drug delivery. Maintains pressure of 13.5 psi.

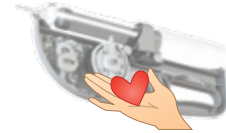

## Advanced technology

All parts are designed to fit & work together within the syringe driver. The Negator motor is the 'heart' of the FREEDOM60<sup>®</sup>.



## Quality assured

Made in the USA with the highest level of quality control for defect-free products. No maintenance needed.




## DynEQ<sup>™</sup> (Dynamic Equilibrium)

If the pressure in your sites increases, the system will immediately respond by automatically slowing down.



## Portable

No batteries or electricity needed. Comes with a travel pouch to promote patient independence.

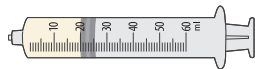

## Standard 60ml capacity

Easy to handle and purge of air. Can be filled from 1ml to 60ml.



## Rugged

Robust design with a molded ABS plastic housing. No sensitive electric parts to get jumbled.



## Spring-loaded syringe holder

Locks the syringe into the driver to maintain proper syringe orientation, prevents inadvertent removal during administration.



## Recessed on/off switch

No unintentional on or off.



## Large wind knob

Ergonomic, easy to grasp & turn, disengages when the FREEDOM60<sup>®</sup> is ON, cannot be over-wound.



## Syringe shield

Easy to check infusion progress, protects the syringe.

**NOTE:** Use ONLY RMS Freedom System accessories with the FREEDOM60<sup>®</sup> Syringe Driver. Use of any other accessory may cause over delivery of fluids, unknown flow rates or additional site complications such as pain, swelling and redness.

## Related Freedom Integrated System Products:



### FreedomEdge<sup>®</sup> Syringe Infusion System

For smaller infusion doses, 20ml & 30ml syringes



### HlgH-Flo Subcutaneous Safety Needle Sets<sup>™</sup>

24G & 26G, 4 - 14mm lengths, 1-6 leg sets



### RMS Precision Flow Rate Tubing<sup>™</sup>

For flow rates from F0.5 - F2400



### 24" Extension Tubing Set (FEXT)

For those who need longer tubing length



rmsmedicalproducts.com 845-469-2042 800-624-9600

The FREEDOM60<sup>®</sup> and FreedomEdge<sup>®</sup> Syringe Infusion Systems, RMS Precision Flow Rate Tubing<sup>™</sup>, HlgH-Flo Subcutaneous Safety Needle Sets<sup>™</sup>, and DynEQ<sup>™</sup> are registered trademarks of RMS Medical Products and are compliant with Medical Device Directive 93/42/EEC. RMS Medical Products is ISO 13485 certified. © 2017 RMS Medical Products; All Rights Reserved. RHL\_INFG\_F60benefits\_03-10005vC